

SAFETY INSTRUCTIONS

DANGER: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

- For adult use only
- Use only on a dry clean surface of the body over clothing by lightly pressing and moving across the body for approximately 60 seconds per region
- Only use the ELEEELS percussive massage device on the body's soft tissue as desired. Do not use on the head or any hard and bone areas of the body. If any pain occurs during use, stop immediately.
- Do not place any objects into the ventilation holes of the ELEEELS percussive massage device.
- Do not allow the device to come into contact with water.
- Do not drop the device.
- Only recharge with the supplied AC adapter.
- Do not tamper or alter the device in any way.
- Never leave the device operating or charging unattended



DISCLAIMER

PLEASE DO NOT USE THE ELEEELS PERCUSSIVE MASSAGE DEVICE, OR ANY OTHER VIBRATION DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM YOUR DOCTOR, SPECIALIST OR PHYSIOTHERAPIST, especially if you are pregnant, have diabetes with complications (such as neuropathy or retinal damage), wear a pacemaker, have had recent surgery, suffer from epilepsy or migraines, herniated disks, spondylolisthesis, spondylo-

sis, or spondylosis, have had recent joint replacements or IUD'S, metal pins or plates, or you have any general concerns about your physical health. Frail Individuals and children should be accompanied by an adult when using any percussive massage devices.

WARRANTY

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN MATERIAL OR WORKMANSHIP, ELEEELS WILL REPAIR OR REPLACE THE PRODUCT, OR NECESSARY COMPONENTS, FREE OF CHARGE

THIS WARRANTY EXCLUDES:

- (a) damage caused by accident, abuse, mishandling, or transport;
- (b) units subjected to unauthorized repair;
- (c) units not used in accordance with ELEEELS care instructions;
- (d) damage exceeding the cost of the product;
- (e) deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) failure to show the dated proof of purchase




Unit 5, 15/F, Global Gateway Tower,
63 Wing Hong Street, Cheung Sha Wan,
Kowloon, Hong Kong
www.eleels.com



X2PRO
Percussive Massage Device

We have developed the ELEEELS percussive massage device to relieve pain and promote circulation, which helps to accelerate recovery and prevents muscle damage.

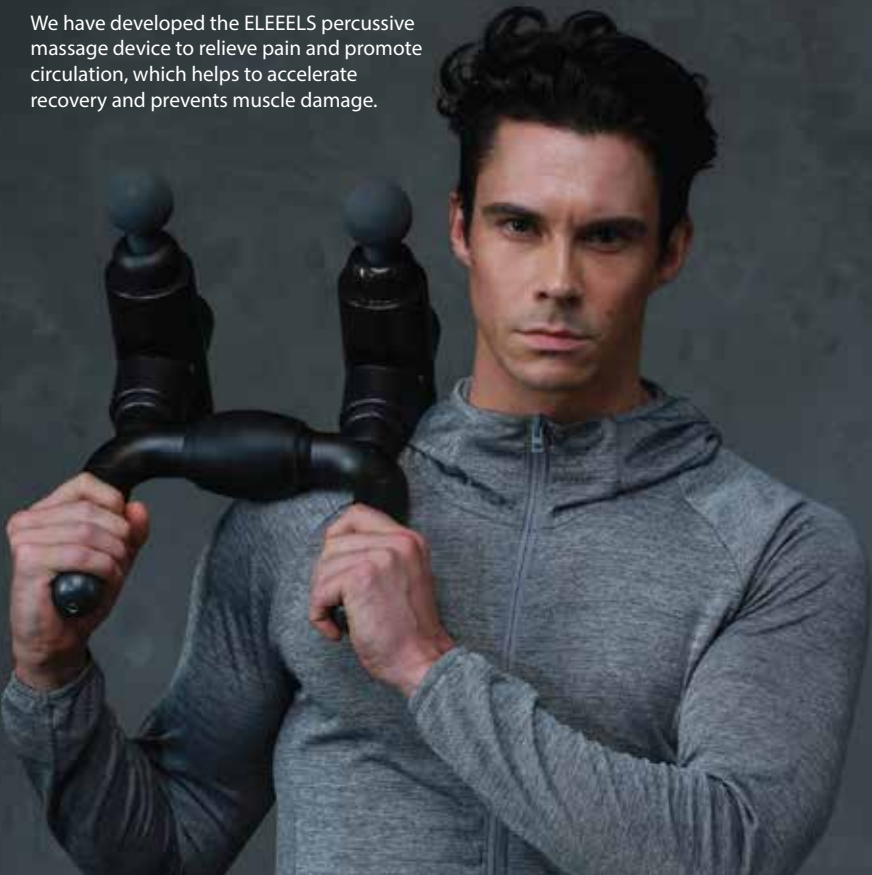

SPECIFICATIONS




1. High Torque and High Speed
2. 6 Speed Light Indicator
3. 6 Speed +/- Controller
4. Quick Start Switch
5. Flexible Ball-joint (Allows free movement of the Twin Heads)
6. Quiet Motion Technology™

1. Powerful Brushless High-Torque Motor
2. 6 Speed Settings for Optimal Pressure
3. ON / OFF Power Switch
4. Interchangeable Applicators (x4)
5. ELEEELS universal adaptor included

Vibration Frequency:
 Comfort Mode: 900rpm, 1200rpm, 1500rpm
 Sports Mode: 1800rpm, 2400rpm, 3200rpm

Weight: 4.4 lbs / 2 kg
 Charge: 100-240V 50/60Hz 4.0A

Interchangeable Applicators

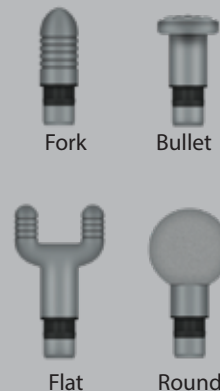

USAGE INSTRUCTIONS

PLEASE READ PRIOR TO USE

- Plug in the supplied 24V 4A AC Adaptor into a mains socket, the LED circle will light up on both handles to indicate power supply.
- For a quick start, turn on the device by touching the +/- speed setting switch on the bottom

right hand side. Pressing the "+" button will enable quick start sports mode and pressing the "-" button will enable quick start comfort mode. Once the desired mode has been selected, simply adjust the speed settings by pressing the +/- buttons on the main switch.

- The massage gun will stop operation after a 10 minute session

Sports Mode level 4-6:

Press the bottom "+" button once to turn device on and into level 4 of sports mode. Pressing the main +/- button can adjust the speed to level 5 and 6 respectively. The blue LED lights will indicate the corresponding power level.

Comfort Mode level 1-3

Press the bottom "-" button once to turn device on and into level 1 of comfort mode. Pressing the main +/- buttons can adjust the speed to level 2 and 3 respectively. The blue LED lights will indicate the corresponding power level.

MAINTENANCE, CLEANING & STORAGE

- Use a slightly damp towel to clean the device and a soft cloth to dry it.
- Do not clean the device with chemicals or alcohol.
- Only use the ELEEELS supplied AC adaptor.

WARNINGS & CAUTION

- Do not immerse in water. Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than one hour. Allow device to rest for 10 minutes before re-using.
- For adult use only. Do not use if injured. Consult your doctor before using this product.