EN

ELEEELS

# A1

Cordless Air  
Compression Leg  
Massage Device

Instruction Manual

**CORDLESS/RECHARGEABLE**  
2000mAh Lithium-ion Battery

**LIGHTWEIGHT DESIGN**  
1.1 lbs. / 500 g per piece

We have developed the ELEEELS cordless air compression leg massage device to massage calf muscle to relieve fatigue and improves circulation, grips and kneads your ankles and heels to push circulation flow upwards through your whole calf, provides rhythmic air compression to strongly squeeze your feet for ultimate relaxation, and grips your shins to relieve pain.

## TECHNOLOGY SPECIFICATIONS

All-in-one Smart  
Controller

6 massage zones

Effective Air Pressure:  
20 kPa - 60 kPa

3 Individual Air Bags


## 3 Different Massage Modes

Rapid (10 mins)  
Gentle (20 mins)  
Dynamic (15 mins)

Power Indicator High  
Low

Battery Indicator  
/ Massage Mode  
Indicator

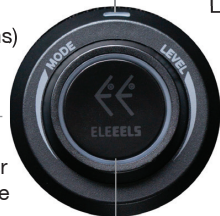- 1) 3 individual air bags
- 2) 6 massage zones
- 3) 2000mAh Lithium-ion Battery  
(up to 3 hours of use per charge)
- 4) Velcro straps
- 5) All-in-one Smart Controller
- 6) 2 power levels
- 7) 3 massage modes for optimal pressure
- 8) Cordless and lightweight (500 g per piece)

Effective Air Pressure:  
Ranges from 20 kPa to 60 kPa

Weight: 1.1 lbs. / 500 g per piece  
Charging Port: USB Type-C

## CHARGING

Fully charge the battery before first use. To charge the battery:


- 1) Connect the USB Type-C cable to the charging port.
- 2) The charging speed depends on the adapter used. It is suggested to use an adapter with at least 5V/2A to fully charge both massage devices in approximately 3 hours.
- 3) During charging, the LED light battery indicator flashes in green slowly. When the green light stops flashing, the device is fully charged.
- 4) The LED light battery indicator shows the battery level of A1:
  - Red: 0-10%
  - Green: 10-100%

## PREPARATION (SLEEVE)


- 1) Please choose the right sleeve first

- 2) Wraps the bands around the heel

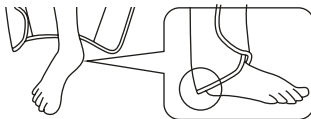- 3) Aligns calf with the upper part of A1

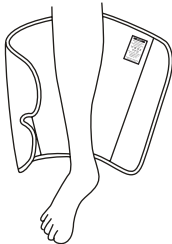
- 4) Wraps at the ankle first, and then the calf. Uses the Velcro straps to adjust the tightness.


- 5) Follows the same procedures for the left sleeve

**PREPARATION (PLANTAR)**

1) Please choose the right sleeve first



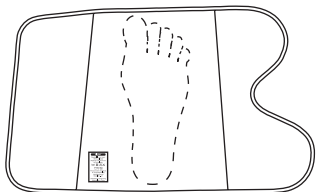
2) Aligns heel with the lower part of A1

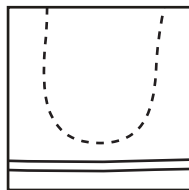3) Wraps at the instep first, and then at the toes. Uses the Velcro straps to adjust the tightness

4) Follows the same procedures for the left sleeve. Please do not walk during plantar massage treatment

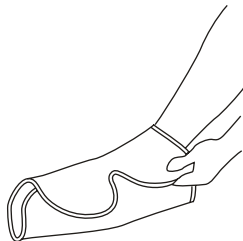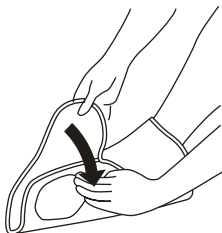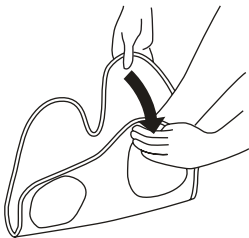**RoHS CE FC**

## OPERATING

- First, presses the all-in-one smart controller twice to switch on.
- Next, switches between power level of “Low” and “High” by rotating the all-in-one smart controller counter clockwise. The “Low” level is indicated by a dimmed power indicator, while the “High” level is indicated by a brightened power indicator.
- Then, switches between “Rapid”, “Gentle”, and “Dynamic” by rotating the all-in-one smart controller clockwise.
- “Rapid” is a 10-minute massage treatment focusing on massaging calf muscle for rapid fatigue relief and improved circulation. The blue LED light massage mode indicator turns on.
- “Gentle” is a 20-minute massage treatment focusing on relaxing the muscle and body comfortably. The purple LED light massage mode indicator turns on.
- “Dynamic” is a 15-minute massage treatment focusing on calf muscle and plantar massage, and relaxation in the same treatment. The white LED light massage mode indicator turns on.
- When the battery level is lower than 10%, the LED light battery indicator turns to red. Operation continues for a few minutes, and the device turns off.
- The device shuts off automatically after each treatment. Alternatively, presses the all-in-one smart controller twice to turn off the device.

## MAINTENANCE, CLEANING & STORAGE

- Use a slightly damp towel to clean the device and a soft cloth to dry it.
- To clean the inner washable layer, detach from Velcro straps and wash at 30°C / 86°F. Do not use bleach. Do not tumble dry.
- Use the ELEEEELS supplied USB Type-C cable.
- Store ELEEEELS A1 in a clean and dry location.

## WARNINGS & CAUTION

- Do not operate continuously for more than one hour. Allow device to rest for 10 minutes before re-using.
- For adult use only. Do not use if injured. Consult your doctor before using this product.

## **SAFETY INSTRUCTIONS**

**DANGER: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:**

- For adult use only.
- If any pain occurs during use, stop immediately.
- Do not allow the device to come into contact with water.
- Do not drop the device.
- Recharges with the USB Type-C cable.
- Do not tamper or alter the device in any way.
- Never leave the device operating or charging unattended.

## **DISCLAIMER**

Please do not use the ELEEELS A1 without first obtaining approval from your doctor if any of the following apply: recent joint replacements, metal pins or plates or you have any concerns about your physical health. Do not use A1 if you feel pain, discomfort, or have a wound at or near the site of application.

## **WARRANTY**

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN MATERIAL OR WORKMANSHIP, ELEEELS WILL REPAIR OR REPLACE THE PRODUCT, OR NECESSARY COMPONENTS, FREE OF CHARGE

## **THIS WARRANTY EXCLUDES:**

- (a) damage caused by accident, abuse, mishandling, or transport;
- (b) units subjected to unauthorized repair;
- (c) units not used in accordance with ELEEELS care instructions;
- (d) damage exceeding the cost of the product;
- (e) deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) failure to show the dated proof of purchase

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