# R2 Lumbar Traction Device

## Instruction Manual

### BACK CARE SOLUTION

Dynamic Stretching  
Vibration Massage  
Thermal Therapy  
Electrical Muscle Stimulation


### OPERATING

- 1) Press auto mode button to turn on the device. Press and hold the button again to turn off.
- 2) Press the auto mode button to switch between four different auto modes: 01: vibration + air inflation, 02: heat + vibration + air inflation, 03: automatic air inflation, 04: heat + air inflation
- 3) Press stretch control button, vibration control button, or heating control button to turn on three different functions separately. The auto mode will stop.
- 4) Press stretch control button, the waist will be lifted automatically to start dynamic stretching. Press the button again once comfortable height is reached. An alarm will sound once you have reached maximum lift. Press and hold the button to release the pressure. The rating bar will blink when releasing and stop blinking when it reaches the original position.
- 5) Press vibration control button to turn on the vibration massage, press and hold the button to turn off. Press the vibration control button to switch the strength between low, medium, and high.
- 6) Press the heating control button to turn on the thermal therapy, press and hold the button to turn off. Press the heating control button to switch the heat between low, medium, and high. The rating bar will blink when reaching the maximum temperature.
- 7) Press the electrical muscle stimulator button to turn on EMS, press and hold the button to turn off. Press the electrical muscle stimulator button to switch between different types of EMS massage. Press the EMS plus button to increase the massage strength, and press the EMS minus button to decrease. The level of the strength will be presented in the touch screen controller. 20 is the highest level.
- 8) To insert an external connector into the electrical pulse output, attach the pads to the corresponding location you would like to massage (arms, legs etc.).
- 9) Reattach pulse pads to protective paper after used.

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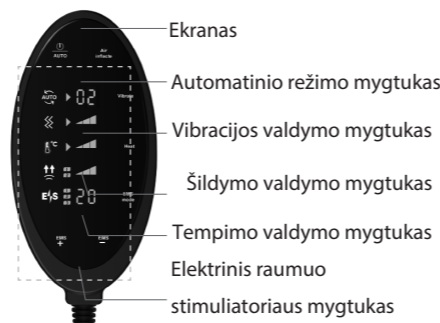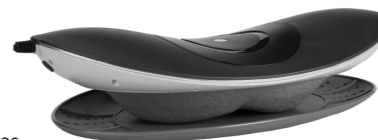
1. The base plate must be placed beneath the device
2. The plate was designed to improve the overall experience when using the traction device on a soft surface
3. You can operate the traction device with or without the plate when the device is placed on hard surface
4. Do not tread on the plate or use it for other purposes

The traction strength will be moderate when the device is laid on this side

The traction strength will be stronger when the device is laid on this side

### USE OF SILICON ELECTRODE SLICES


- Insert one end of the 2-core connector wires into the electrode port and connect the other end to the electrode slices.
- The electrode slices should be clean on their surfaces without any dust, oily or sticky substance. If their surfaces become dirty, flush with clean water and then air dry them before repeated use.
- The electrode slices can be used for around 50 times of 20 minutes session. If their viscosity cannot be recovered after several times of cleaning, please replace new ones.
- Only use the electrode slices provided with the device.
- Please dispose electrode slices according to your local environmental regulations rather than discarding them at will.

NOTE: It is recommended to use the EMS massager once or twice a day. Please give your skin a break from the electrode pads in every 20 minutes. Leaving pads on the same skin location may cause skin irritation.

### MAINTENANCE, CLEANING & STORAGE

- Store this device in a moisture-free, dry place when not in use.
- Do not wrap the power cord around the device.
- Do not use the power cord as a handle.
- Keep away from all solvents and harsh detergents.
- Unplug the device and allow it to cool before cleaning. Wipe with a soft, dry cloth as needed.
- Do not submerge the device in water or other liquids.
- Do not use abrasive cleaners, brushes, gasoline, kerosene, glass/furniture polish, or paint thinner to clean the massager.

### PREPARATION

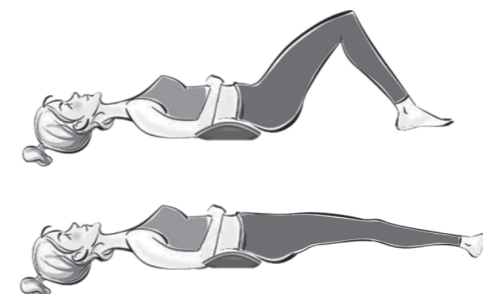- Please put the device under your waist between lumbar no. L1 and L5.
- Please place device on a hard surface or use a hard cushion under the device to keep level for operation.

### FEATURES

- Air traction system with ergonomically designed to fit body curves
- Infrared light and heat therapy with adjustable temperature level
- Vibration massage with two high-power 3000 rpm motors
- Includes one electromagnetic chip and EMS pads for massaging other body parts
- ABS & PC materials supporting more than 200 kg of weight

### NOTE

- The default time for infrared heat therapy is 15 minutes. After a 10-minute break, press the on/off button on the control panel again if you need more cycles.
- In different mode, some buttons for different functions will be locked.
- This device is only used for waist.
- Do not contact skin directly with this device.

### SAFETY INSTRUCTIONS

- 1) Only use the adapter that was supplied with the device.
- 2) Do not let the device fall from high place, and do not keep it around sharp objects.
- 3) Do not use this product when you are sleeping.
- 4) Do not use other medical electronic devices at the same time, such as pacemakers, heart-lung machines or electrocardiographs. If the patient is simultaneously using high frequency surgical equipment or other related equipment, the electrical pulse patches in the device may cause burns or damage the equipment. If the equipment is used near a microwave, the output of the device may become unstable.
- 5) Keep away from children.
- 6) When using the electric pulse function, do not move or change the position.

### DISCLAIMER

Please do not use the ELEEEELS R2 without first obtaining approval from your doctor if any of the following apply: pregnant, have skin sensitivities, heart disease, blood pressure abnormalities, malignant tumors, cerebrovascular diseases, diabetes, recent joint replacements, metal pins or plates or you have any concerns about your physical health. Do not use R1 if you feel pain, discomfort, or have a wound at or near the site of application.

### WARRANTY

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN MATERIAL OR WORKMANSHIP, ELEEEELS WILL REPAIR OR REPLACE THE PRODUCT, OR NECESSARY COMPONENTS, FREE OF CHARGE

### THIS WARRANTY EXCLUDES:

- (a) damage caused by accident, abuse, mishandling, or transport;
- (b) units subjected to unauthorized repair;
- (c) units not used in accordance with ELEEEELS care instructions;
- (d) damage exceeding the cost of the product;
- (e) deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) failure to show the dated proof of purchase

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