ELEEEELS

P1


Piculet
Percussive
Massager

Instruction Manual

We have developed the ELEEEELS percussive massager to accelerate recovery, prevents muscle damage, relieves pain, promotes circulation, and everyday relaxation anytime, anywhere.

P1 at a Glance

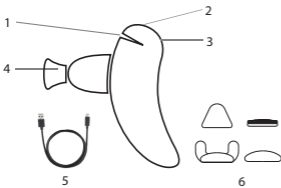- 1) Speed / battery level indicator
- 2) Power / speed control button
- 3) USB-C charging port cover
- 4) Prime CNC aluminium alloy applicator
- 5) USB-C cable
- 6) 4 add-on silicon rubber attachments

Vibration Frequency:

Level 1 (20 strokes per second) / Level 2 (30 strokes per second)
Level 3 (40 strokes per second) / Level 4 (50 strokes per second)

Specifications

- Rated Voltage: 7.4V
- Rated Power: 15.16W
- Input: USB-C 5V/2A
- Battery Capacity: 1000 mAh
- Charging Port: USB-C
- Weight: 277 grams
- Operating Time: around 2.5 hours of use per charge under maximum speed

INSTRUCTIONS

PLEASE READ PRIOR TO USE

CHARGING

- 1) Flip over the charging port cover gently to unveil the USB-C charging port. Connect USB-C cable to the charging port.
- 2) The charging speed depends on what adapter is used. It is suggested to use an USB adapter supporting 5V/2A to fully charge in approximately 2.5 hours.
- 3) During charging, the battery level indicator flashes in green slowly. When the device is fully charged, the battery level indicator turns off.

OPERATING

- First, press the power button for 1.5 seconds to switch on to level 1.
- Then, press the speed control button repeatedly to select the required speed level. The speed indicator shows the speed level:
- Level 1: Pink / Level 2: Purple / Level 3: Blue / Level 4: White
- Massage the desired body part at the desired speed by applying pressure as needed.
- Press the power button for 1.5 seconds to switch off.
- When the massage gun stops operation after a 10-minute session, simply press the power button to start a new one.
- When the battery level is lower than 10%, the battery level indicator flashes in red. Operation stops. Please recharge the battery.

THERMAL PERCUSSIVE MASSAGE - TUTORIAL


Press the power button for 1.5 seconds to switch on. Heat is generated directly from P1's high-speed operation. After 5 minutes of operation (assume working at room temperature), the prime CNC aluminium alloy applicator warms up to around 37°C - 45°C automatically without extra battery or heating plate. CNC aluminum alloy, an ideal conductor, helps to transfer the heat from P1 to the massage area.

The warmup speed of the applicator depends on the speed level. Higher speed level generates more heat energy, which will lead to higher temperatures of the applicator. Enjoy thermal percussive massage by massaging the desired body part.

PRIME APPLICATOR / ADD-ON ATTACHMENTS

You may use the prime CNC aluminium alloy applicator only to enjoy massage on all body parts. For additional massage experiences, please insert an add-on silicon rubber attachment to the prime applicator by align it with the tip of the applicator and push to connect. To remove the attachment, place 1 finger on each side of the attachment and pull it off from the applicator.

Massage different body parts with the prime applicator alone, or with the prime applicator and the corresponding add-on attachment together to fit a wide range of therapeutic and relaxation needs.

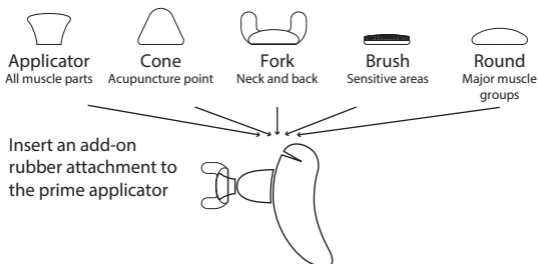MAINTENANCE, CLEANING & STORAGE

- Use a slightly damp towel to clean the device and a soft cloth to dry it.
- Do not clean the device with chemicals or alcohol.
- Only use the USB-C cable supplied by the ELEEEELS for recharging.

WARNINGS & CAUTION

- Do not immerse in water. Keep away from liquids or heat sources.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than one hour. Allow device to rest for 10 minutes before re-using.
- For adult use only. Do not use if injured. Consult your doctor before using this product.

SAFETY INSTRUCTIONS

DANGER: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

- Apply to a clean, dry clothed area of the body by lightly pressing and moving across the body for approximately 60 seconds per area.
- Only use the device on the body's soft tissue as desired. Do not use on head or any hard, bony area of the body. Discontinue use immediately if experiencing any pain.
- Do not drop the device, and do not tamper or alter the device in any way.
- Never leave the device unattended during operating or charging.

DISCLAIMER

Please do not use the ELEEEELS P1 without first obtaining approval from your doctor if you have any of the following conditions: pregnancy, neuropathy, retinal damage, pacemaker user, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements, surgical pins or plates, or you have any concerns about your physical health.

WARRANTY

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN THE MATERIAL OR WORKMANSHIP, ELEEEELS WILL REPAIR OR REPLACE THE PRODUCT; OR FAILED COMPONENT(S), FREE OF

CHARGE. THIS WARRANTY EXCLUDES:

- (a) Damage caused by accident, abuse, or mishandling during transport;
- (b) Damage arising from use of units from unauthorized repair;
- (c) Malfunction of units that are not used in accordance with ELEEEELS care instructions;
- (d) Damage for which the cost of repair exceeds the cost of the product;
- (e) Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) Any product for which a dated proof of purchase cannot be provided.

Please visit below link for instruction manual in other languages:
<http://www.eleeeels.com/im>

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