# X1T

Percussive  
Massage  
Device

## Instruction Manual

### CORDLESS/RECHARGEABLE

24 V Lithium-ion Battery

### LIGHTWEIGHT DESIGN

2.5lbs / 1.1kg

- Red: 0-25%
- Blue: 25-75%
- Green: 75- 100%

### HOW DO I ADJUST THE HANDLE'S ANGLE?

Twist the battery (handle) 180 degrees. A 'click' can be heard when this is completed

### MAINTENANCE, CLEANING & STORAGE

- Use a slightly damp towel to clean the device and a soft cloth to dry it.
- Do not clean the device with chemicals or alcohol.
- Only use the ELEEELS supplied AC adaptor.

### WARNINGS & CAUTION

- Do not immerse in water. Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than one hour. Allow device to rest for 10 minutes before re-using.
- For adult use only. Do not use if injured. Consult your doctor before using this product.


### SAFETY INSTRUCTIONS

DANGER: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

- For adult use only
- Use only on a dry clean surface of the body over clothing by lightly pressing and moving across the body for approximately 60 seconds per region
- Only use the ELEEELS percussive massage device on the body's soft tissue as desired. Do not use on the head or any hard and bone areas of the body. If any pain occurs during

We have developed the ELEEELS percussive massage device to relieve pain and promote circulation, which helps to accelerate recovery and prevents muscle damage.

### SPECIFICATIONS

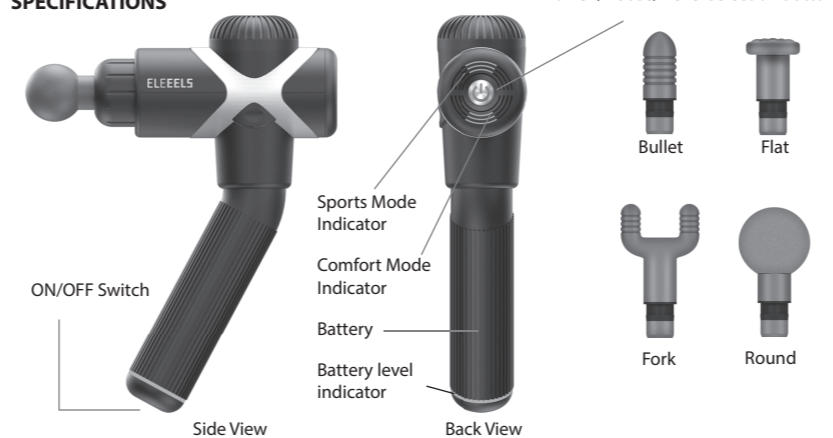- 1 Powerful Brushless High-Torque Motor
- 2 Rechargeable Lithium-Ion Battery (3+ Hrs of Use Per Charge)
- 3 6 Speed Settings for Optimal Pressure
- 4 Cordless and Lightweight (2.5 lbs)
- 5 ON / OFF Power Switch
- 6 Interchangeable Applicators (x4)
- 7 ELEEELS universal adaptor included

Vibration Frequency:  
Comfort Mode: 900rpm, 1200rpm, 1500rpm  
Sports Mode: 1800rpm, 2400rpm, 3200rpm  
Weight: 2.5 lbs / 1.1 kg  
Charge: 110-240V 50/60Hz 28V, DC-1A

use, stop immediately.

- Do not place any objects into the ventilation holes of the ELEEELS percussive massage device.
- Do not allow the device to come into contact with water.
- Do not drop the device.
- Only recharge with the supplied AC adaptor.
- Do not tamper or alter the device in any way.
- Never leave the device operating or charging unattended

### DISCLAIMER

PLEASE DO NOT USE THE ELEEELS PERCUSSIVE MASSAGE DEVICE, OR ANY OTHER VIBRATION DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM YOUR DOCTOR, SPECIALIST OR PHYSIOTHERAPIST, especially if you are pregnant, have diabetes with complications (such as neuropathy or retinal damage), wear a pacemaker, have had recent surgery, suffer from epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylolysis, have had recent joint replacements or IUD'S, metal pins or plates, or you have any general concerns about your physical health. Frail Individuals and children should be accompanied by an adult when using any percussive massage devices.

### WARRANTY

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN MATERIAL OR WORKMANSHIP, ELEEELS WILL REPAIR OR REPLACE THE PRODUCT, OR NECESSARY COMPONENTS, FREE OF CHARGE

THIS WARRANTY EXCLUDES:

### USAGE INSTRUCTIONS

PLEASE READ PRIOR TO USE

### CHARGING

Fully charge the battery for up to 6 hours before first use.

To charge the battery:

- 1) Turn off the main power switch located at the bottom of the handle and connect to the supplied adaptor.
- 2) A green light will flash at the bottom of the handle when charging has begun
- 3) When the green light has stopped flashing, the device has been fully charged

### OPERATING

- Press the main power button at the bottom of the handle to "ON".
- Press the power button located at the rear of the device.
- The massage gun will stop operation after a 10 minute session, simply press the power button to start a new one.

Sports Mode level 1-3:

Press the bottom power button once to turn device on and into level 1 of sports mode. Pressing the main power button can adjust the speed to level 2 and 3 respectively. The blue LED lights will indicate the corresponding power level.

Comfort Mode level 1-3

Hold the main power button for 2 seconds and release to turn device on and into level 1 of comfort mode. Further pressing the main power button can adjust the speed to level 2 and 3 respectively. The blue LED lights will indicate the corresponding power level.

- Massage the desired body part at the desired speed by applying pressure as needed
- The circular light indicator on the battery:

- (a) damage caused by accident, abuse, mishandling, or transport;
- (b) units subjected to unauthorized repair;
- (c) units not used in accordance with ELEEELS care instructions;
- (d) damage exceeding the cost of the product;
- (e) deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) failure to show the dated proof of purchase

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