ELEEEELS

X3

Percussive
Massage
Device

Instruction Manual


CORDLESS/RECHARGEABLE
14.8V Lithium-ion Battery

LIGHTWEIGHT DESIGN
2 lbs. / 0.9 kg

We have developed the ELEEEELS percussive massage device to accelerate recovery and prevents muscle damage so as to relieve pain and promote circulation.

EN

SPECIFICATIONS

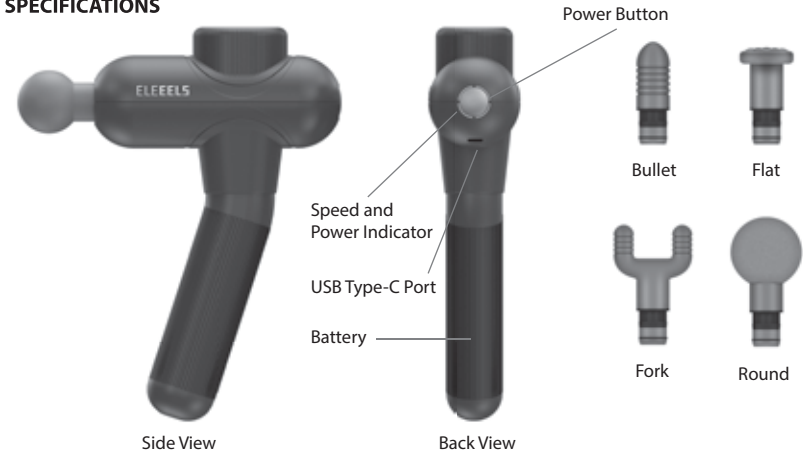- | | |
|--|---|
| 1) High-speed brushless motor (max 3200 rpm) | 9) Interchangeable applicators (x4) |
| 2) Endurance Motion Technology™ | 10) ELEEEELS USB Type-C cable included |
| 3) Rechargeable lithium-Ion battery (up to 10 hours of use per charge) | 11) ELEEEELS carrying case included |
| 4) Quiet Motion Technology™ | Vibration Frequency:
Level 1 (1200 rpm), Level 2 (1800 rpm),
Level 3 (2400 rpm), Level 4 (3200 rpm) |
| 5) 4 speed settings for optimal pressure | Weight: 2 lbs. / 0.9 kg |
| 6) Cordless and lightweight (0.9 kg) | Charging Port: USB Type-C |
| 7) Power button | |
| 8) Speed and power indicator | |

INSTRUCTIONS

PLEASE READ PRIOR TO USE

CHARGING

Fully charge the battery for up to 6 hours before first use. To charge the battery:

- 1) Connect USB Type-C cable to the charging port
- 2) The charging speed depends on what adapter is used. It is suggested to use a adapter supporting Type C (5V/3A), QC 3.0, or PD (9V/2A) to fully charge in approximately 3 hours.

EN

3) During charging, the green LED light power indicator flashes slowly. When the green light stops flashing, the device is fully charged.

4) The green LED light power indicator shows the battery level:

- Red: 0-10%
- Green (level 1): 11- 25%
- Green (level 2): 26-50%
- Green (level 3): 51-75%
- Green (level 4): 76-100%

OPERATING

- First, press the power button for 2 seconds to switch on. The green LED light power indicator will light to show the battery level.
- Then, press the power button repeatedly to select the required level. The currently selected power level is shown by the blue LED light operating indicator.
- When the battery level is lower than 10%, the power indicator turns red. Operation continues, but it is recommended to recharge the battery.
- When the device stops operation for 2 minutes, the green LED light power indicator will flash for 2 times before the device shuts off automatically.
- Massage the desired body part at the desired speed by applying pressure as needed.
- When the massage gun stops operation after a 10-minute session, simply press the power button to start a new one.

MAINTENANCE, CLEANING & STORAGE

- Use a slightly damp towel to clean the device and a soft cloth to dry it.
- Do not clean the device with chemicals or alcohol.
- Only use the USB Type-C cable supplied by the ELEEEELS for recharging.

WARNINGS & CAUTION

- Do not immerse in water. Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than one hour. Allow device to rest for 10 minutes before re-using.
- For adult use only. Do not use if injured. Consult your doctor before using this product.

SAFETY INSTRUCTIONS

DANGER: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

- For adult use only.

EN

- Apply to a clean, dry clothed area of the body by lightly pressing and moving across the body for approximately 60 seconds per area.
- Only use the ELEEEELS percussive massage device on the body's soft tissue as desired. Do not use on the head or any hard, bony areas of the body. Discontinue use immediately if experiencing any pain.
- Do not allow the device to come into contact with water.
- Do not drop the device.
- Only recharge with the supplied USB Type-C cable.
- Do not tamper or alter the device in any way.
- Never leave the device operating or charging unattended .

DISCLAIMER

Please do not use the ELEEEELS X3, or any massage device without first obtaining approval from your doctor if you have any of the following conditions: pregnancy, neuropathy, retinal damage, pacemaker user, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements, surgical pins or plates, or you have any concerns about your physical health.

WARRANTY

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN THE MATERIAL OR WORKMANSHIP, ELEEEELS WILL REPAIR OR REPLACE THE PRODUCT; OR FAILED COMPONENT(S), FREE OF CHARGE

THIS WARRANTY EXCLUDES:

- (a) Damage caused by accident, abuse, or mishandling during transport;
- (b) Damage arising from use of units from unauthorized repair;
- (c) Malfunction of units that are not used in accordance with ELEEEELS care instructions;
- (d) Damage for which the cost of repair exceeds the cost of the product;
- (e) Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) Any product for which a dated proof of purchase cannot be provided.

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