ELEEEELS

**P2** | Pocket  
Massage  
Gun

Instruction Manual


**CORDLESS/RECHARGEABLE**  
7.4 V Lithium-ion Battery

**LIGHTWEIGHT DESIGN**  
0.79 lbs. / 360 g

We have developed the ELEEEELS massage gun to accelerate recovery, prevents muscle damage, relieves pain, promotes circulation, and everyday relaxation anytime, anywhere.

EN

### SPECIFICATIONS

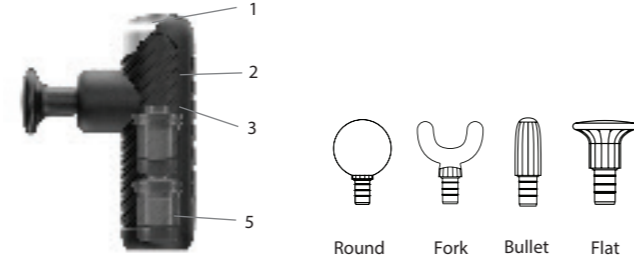- 1) Durable and lightweight with aluminium alloy body inside (360 grams)
- 2) High-torque and high-speed 38 mm brushless motor (delivers up to 50 strokes per second)
- 3) Rubberized surface for a no-slip grip
- 4) FOUR speed settings and FOUR applicators to fit a wide range of therapeutic and relaxation needs
- 5) Rechargeable lithium-ion battery with 7.4V output voltage (up to 6 hours of use per charge)
- 6) Mini Piston Technology

Vibration Frequency:  
Level 1 (30 strokes per second), Level 2 (37 strokes per second)  
Level 3 (43 strokes per second), Level 4 (50 strokes per second)

Weight: 0.79 lbs. / 360 g  
Charging Port: USB Type-C


### INSTRUCTIONS

PLEASE READ PRIOR TO USE

### CHARGING

- 1) Connect USB Type-C cable to the charging port
- 2) The charging speed depends on what adapter is used. It is suggested to use an adapter supporting Type C (5V/3A), QC 3.0, or PD (9V/2A) to fully charge in approximately 3 hours.

EN

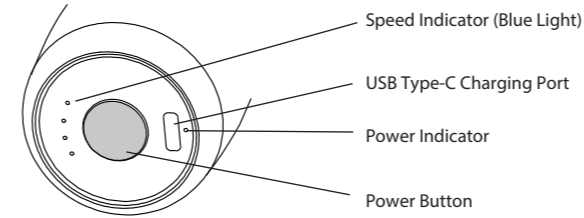- 3) During charging, the power indicator is in red solid light. When the green solid light is on, the device is fully charged.
- 4) The power indicator shows the battery level:
  - Red (flashes 1 time): 0-10%
  - Green (flashes 1 time): 11- 30%
  - Green (flashes 2 times): 31-70%
  - Green (flashes 3 times): 71-99%
  - Green (solid light): 100%

### OPERATING

- First, press the power button for 1.5 seconds to switch on. The power indicator will light to show the battery level.
- Then, press the power button repeatedly to select the required speed level. The currently selected speed level is shown by the blue speed indicator.
- Press the power button for 1.5 seconds to switch off.
- Massage the desired body part at the desired speed by applying pressure as needed.
- When the massage gun stops operation after a 10-minute session, simply press the power button to start a new one.

### MAINTENANCE, CLEANING & STORAGE

- Use a slightly damp towel to clean the device and a soft cloth to dry it.
- Do not clean the device with chemicals or alcohol.
- Only use the USB Type-C cable supplied by the ELEEEELS for recharging.

### WARNINGS & CAUTION

- Do not immerse in water. Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than one hour. Allow device to rest for 10 minutes before re-using.
- For adult use only. Do not use if injured. Consult your doctor before using this product.

EN

### SAFETY INSTRUCTIONS

**DANGER: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:**

- Apply to a clean, dry clothed area of the body by lightly pressing and moving across the body for approximately 60 seconds per area.
- Only use the device on the body's soft tissue as desired. Do not use on the head or any hard, bony areas of the body. Discontinue use immediately if experiencing any pain.
- Do not drop the device.
- Do not tamper or alter the device in any way.
- Never leave the device unattended during operating or charging.

### DISCLAIMER

Please do not use the ELEEEELS P2, or any massage device without first obtaining approval from your doctor if you have any of the following conditions: pregnancy, neuropathy, retinal damage, pacemaker user, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements, surgical pins or plates, or you have any concerns about your physical health.

### WARRANTY

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN THE MATERIAL OR WORKMANSHIP, ELEEEELS WILL REPAIR OR REPLACE THE PRODUCT; OR FAILED COMPONENT(S), FREE OF CHARGE

THIS WARRANTY EXCLUDES:

- (a) Damage caused by accident, abuse, or mishandling during transport;
- (b) Damage arising from use of units from unauthorized repair;
- (c) Malfunction of units that are not used in accordance with ELEEEELS care instructions;
- (d) Damage for which the cost of repair exceeds the cost of the product;
- (e) Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) Any product for which a dated proof of purchase cannot be provided.

Please visit below link for instruction manual in other languages:  
<http://www.eleeeels.com/im>

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