S1

Revival Hot
Stone Spa
Collection

ELEEELS

Instruction Manual

EN

Find your personal gateway within the comfort of your own home with the S1 Revival Hot Stone Spa Collection.

ELEEELS S1 at a Glance

- 1) 8 hot stones made from basalt for heat therapy and hot stone massage
- 2) Ambient light therapy control
- 3) ON/OFF switch for temperature control
- 4) Charging port
- 5) Aromatherapy massage oil

Weight: 4.41 lbs. / 2 kg
 Input: 100-240V 50/60Hz 1.0A
 Max Output: 24.0V 2.0A 48.0W


3 Hot Stone Sizes

Sacrum Stone

Diameter: 3.25 inches
 Designed for treatments on lower back (sacrum), abdomen (belly), or between shoulders

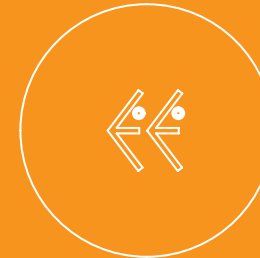Standard Stone

Diameter: 2.5 inches
 Designed for treatments on arms, thighs, calves, palm of hands, arch of feet and back of the neck.

Facial Stone

Diameter: 2 inches
 Designed for treatments on face, hands, arms, and neck.

OPERATING - HEAT THERAPY AND HOT STONE MASSAGE

- Switch the power button up to turn on the console and ambient light therapy with blinking green lights. There are 2 heat settings. Switch the power button to the middle for low heat, and switch the power button to the top for high heat. The temperature for low heat is about 48°C, while the temperature for high heat is about 55°C.

- Always start with low heat for first-time use.

- Place the hot stones on the console to heat up. Hot stone massage will be ready in about 10 minutes. For higher temperature, heat the hot stones for a longer time.

- **Always touch the stones first to feel the temperature before holding it firmly on hands, and applying on body.**


- Place the hot stones either on nape, back, waist, arms, shoulders, or legs of your partner to enjoy heat therapy.
- To enjoy hot stone massage, apply a small amount of ELEEELS Care Massage Oil on the skin of your partner before rubbing hot stones on different parts of the body.
- Pick up any of the hot stone(s) and rub them on target areas of the body.
- Place the hot stones back on the console to re-heat or keep warm if necessary.
- Switch the power button down to turn off the console. Alternatively, the console will turn off automatically after a 90-minute session. Please also switch the power button down as well.


OPERATING - AMBIENT LIGHTS


The ambient light will automatically start once the power is switched on. The default mode of the ambient lights is blinking green lights.

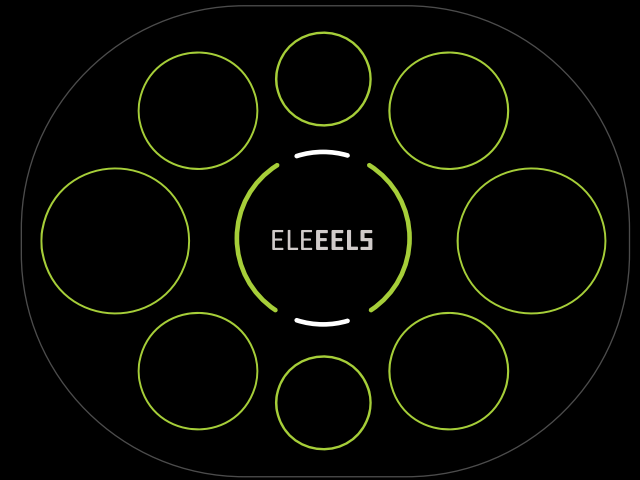
Tap the top button once for slow blinking rainbow lights.

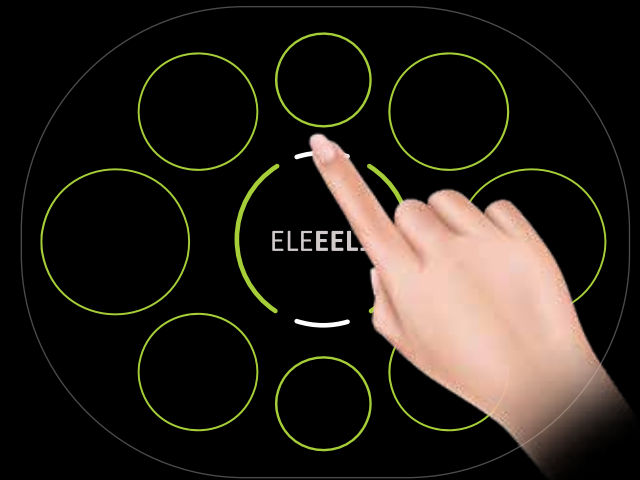
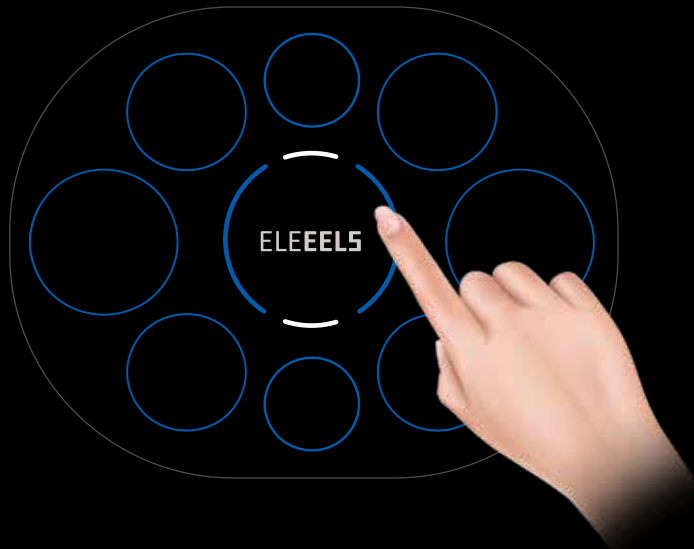
Swipe through the right LED light ring for solid rainbow lights.

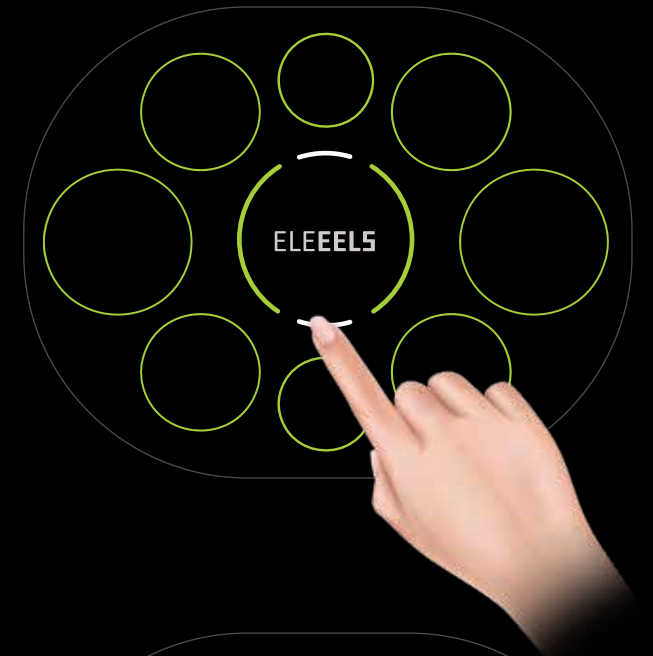Tap the bottom button alternatively to adjust the three-level light intensity (only for solid color lighting).

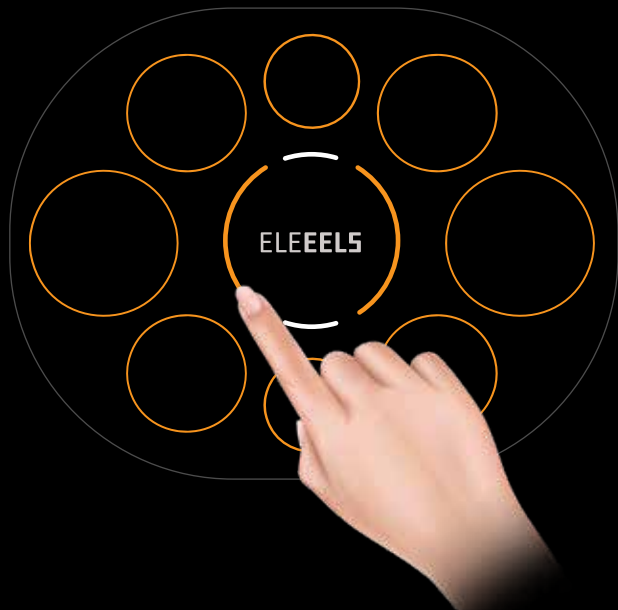
Swipe through the left LED light ring for blinking rainbow lights.

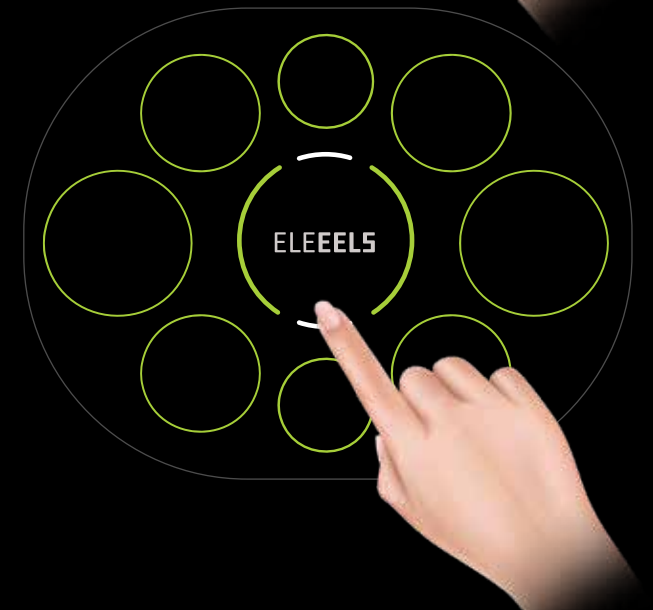Tap and hold the bottom button for 3 seconds to turn on and off the ambient lights.

WARNINGS & CAUTIONS

- This device is specifically intended for body massage using hot basalt stones.
- This device and its components are not designed to be used by children. Keep away from children.
- Gently place the hot stones on the console.
- **Always touch the stones first to feel the temperature before holding it firmly on hands, and applying on body.**
- **Be aware that the heating plate can be extremely hot during operation, and remains hot for a while after the device has been turn off. NEVER touch the heating plate when it is hot.**

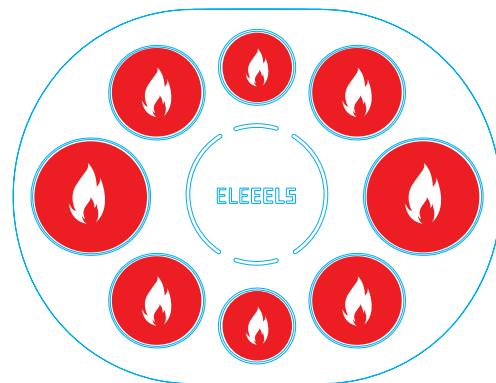WARNINGS & CAUTIONS

- Stop the massage immediately if you feel uncomfortable.
- Do not use the device in humid rooms.
- Always switch off the product and remove the plug from the wall socket when not in use and before cleaning.
- Ensure the product is switched off before disconnecting it from the power supply.
- Avoid touching the plug or power switch with wet or damp hands.
- Do not place the device next to a heat source, water, or other liquids.
- Do not move the appliance when it is in use.
- Do not place any objects or pour liquid on the console.
- If this product stops working, switch off and unplug the product immediately. Seek professional advice before attempting to re-use this product.
- Do not remove screws or attempt to disassemble.
- Do not drop the appliance.
- Do not tamper or alter the appliance in any way.

MAINTENANCE & CLEANING

- Use a soft towel to clean the device. Add a minimal amount of detergent to clean the massage oil.
- Do not clean the device with chemicals or alcohol.
- Only use the AC adapter supplied by ELEEEELS.

DISCLAIMER

Please do not use the ELEEEELS S1 without first obtaining approval from your doctor if you are pregnant, pacemaker user, or diagnosed with any of the following conditions: heart disease, malignant tumor, scoliosis or abnormality, hypertension, osteoporosis, thrombosis, severe aneurysm, dermatitis and other skin infections (including inflammation of subcutaneous tissue), wounds, and severe circulatory disorders.

Our hot stone is a natural product and has been hand polished as such contains slight variations in color, veining, texture, contour, or shading.

WARRANTY

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN THE MATERIAL OR WORKMANSHIP, ELEEEELS WILL REPAIR OR REPLACE THE PRODUCT; OR FAILED COMPONENT(S), FREE OF CHARGE.

THIS WARRANTY EXCLUDES:

- (a) Damage caused by accident, abuse, or mishandling during transport;
- (b) Damage arising from use of units from unauthorized repair;
- (c) Malfunction of units that are not used in accordance with ELEEEELS care instructions;
- (d) Damage for which the cost of repair exceeds the cost of the product;
- (e) Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) Any product for which a dated proof of purchase cannot be provided.

Please visit below link for instruction manual in other languages:

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